

**Cimbria Heid GmbH**  
Heid Werkstrasse 4  
A-2000 Stockerau, Austria  
Telefon: (43) -2266-6990  
Telefax: (43) -2266-65590  
E-mail : [heid@cimbria.at](mailto:heid@cimbria.at) or [wbo@cimbria.at](mailto:wbo@cimbria.at)  
Internet: [www.cimbria.com](http://www.cimbria.com) or [www.cimbriakaack.com](http://www.cimbriakaack.com)




## Sample Screens - Sieving Software and Hardware

### Original Sample Screens from Kaack with wooden frames



No. 10 = Ø 3,97 mm  
No. 11 = Ø 4,37 mm  
No. 12 = Ø 4,76 mm  
No. 13 = Ø 5,16 mm  
No. 14 = Ø 5,56 mm  
No. 15 = Ø 5,95 mm  
  
No. 16 = Ø 6,35 mm  
No. 17 = Ø 6,75 mm  
No. 18 = Ø 7,14 mm  
No. 19 = Ø 7,54 mm  
No. 20 = Ø 7,94 mm  
Blind bottom and Cover

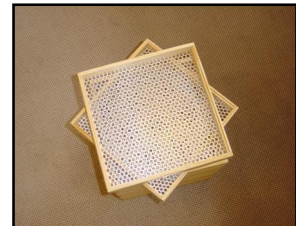

### Software - Sample Check V2.00

- Master Data
- Order Number – LOT/ID-Number – Region – Quality defaults
- Screen weight /weight acquisition
- Sample comparison / analysis
- Data logging
- Link to contract administration – Network Management

### Laboratory Vibrating sample screen machine LS 500

With linear motion production an intensive shaking of sample screens with dimensions of 300mm x 300mm. With compact housing and 4 rubber feet. Dimensions of the shaking table 540 mm x 380 mm with mounted holding device with quick action levers for sample screen.

Driven by speed stabile A/C motor with V-belt drive and fixed transmission ratio. Built-in mass equalizer for a stable stand and guaranteeing an exact shaking. Silent running by means of lifetime

greased bearings for the driving elements as well as maintenance free plastic bearings for the moving parts. With built-in electronic timer allow to pre-set the wanted operation time with high repeating exactness.

